# RENPHO

SMART HEALTHY LIVING

Model:RF-EM001R

## Eye Massager

**US Importer: JOICOM CORPORATION**

14129 The Merge Street, Building 3 Unit A, Eastvale, CA 92880

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**Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.**

Qianhai Complex A201, Qianwan Road 1, Qianhai Shenzhen-Hong Kong Cooperation Zone, Shenzhen, P.R.China

**EU REP C&E Connection E-Commerce (DE) GmbH**

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---

EN

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

DE

Children shall not play with the appliance.

FR

Cleaning and user maintenance shall not be made by children without supervision.

IT

ES

Please read and comprehend all the contents of this manual and operate the product according to the instructions.

1. **DO NOT** use if you have had an eye operation, or if you have any eye diseases such as glaucoma, cataract, and amotio retinae.
2. **DO NOT** clean the product with water; clean it with cotton soaked in alcohol
3. This product has a hot surface; people sensitive to heat should take extra care when using this.
4. Juveniles, disabled people, and the elderly should only use this product when properly supervised.
5. The product is intended for personal use; it is not intended to replace medical care and medical devices.
6. **DO NOT** use if you have any swelling in the eye area.
7. **DO NOT** disassemble the product, **DO NOT** impact or squeeze its interior lithium battery or expose it to fire;
8. **DO NOT** use if the battery has been immersed in water
9. **DO NOT** use in high temperatures or humid environments. Avoid direct exposure to the sun.
10. Consult your doctor if you have sensitive skin, allergic dermatitis, or orthopedic bio materials.
11. Remove contact lenses, hairpins, and hair accessories before using.
12. Start the eye massager after putting it on in order to avoid damage to the air bag due to over pressure.
13. Only use the charger (5V 1A, not included) to recharge the eye massager.
14. You must have your eyes closed while massaging.
15. **DO NOT** use sharp objects to scratch this product

## Product Structure

---


EN

DE

FR

IT

ES
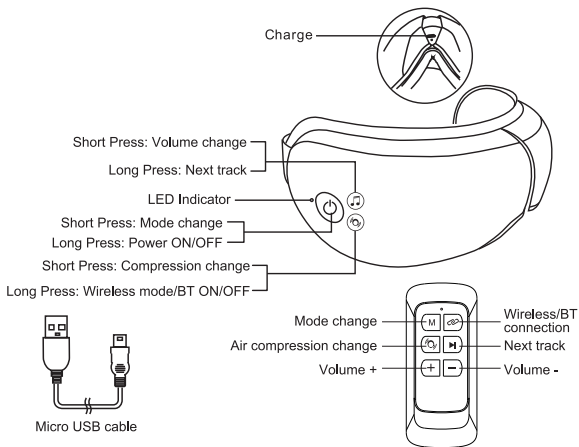
## Technical parameters

---

Product name: Eye Massager

Model: RF-EM001R

Input: 5V  $\equiv$  1A

Rated voltage: 3.7V

Rated power: 5W

Alkaline battery of remote: 23A/12V

## What's included

---

1 x RENPHO Eye Massager

1 x USB Cable

1 x User Manual

1 x Remote Control

1 x Storage Bag

1 x Alkaline Battery

## Wearing the device

Align the cloth cover to your eyes, put on the device, adjust the elastic band, adjust the device to the most comfortable position. Hold the power button for 3-5 seconds to turn on the eye massager

EN

DE

FR

IT

ES

Attention

Remove any glasses, contact lenses, hairpins and hair accessories before use to avoid discomfort. Do not operate the device with wet hands.


## During operation

Default Setting: Air Pressure + Heat + Music

Default Timer: 15 minutes

Press power button on the device, or “M” on the remote to change mode:


Hold power button 3-5 seconds to turn on/off the eye massager. Gently press power button to change modes. Audio voice will confirm mode change.
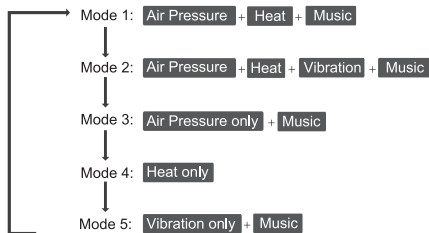
# Charging

---

1. Use the 5V mobile charger(not included) for charging your eye massager.
2. Insert USB cable into the socket of the mobile charger.
3. Softly insert the micro USB plug into the socket of the device, making sure the micro USB is inserted in the right direction.
4. Plug the charger into power.
5. While charging, the red LED indicator will show. When done charging, the LED will be green. The device will not work while charging.
6. The charge duration will be 2-3 hours. Disconnect the power adapter after charging is complete.

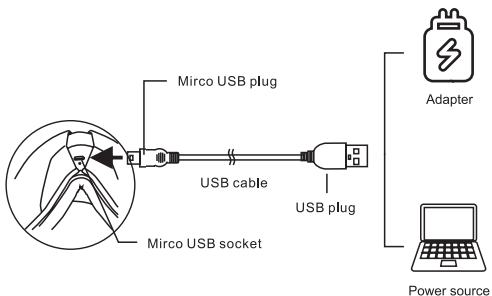Note: The adapter is not included.

EN

DE

FR

IT

ES

EN

## Volume adjustment

Short press " " to adjust the volume among: Off-Volume 1-Volume 2-Volume 3-Volume 4-Volume 5(max)-Off.

DE

## Wireless/BT Pairing

The product automatically enters the Wireless/BT pairing mode when it is turned on, and the blue indicator light flashes slowly.

Turn on the Bluetooth on your phone and found the Bluetooth name "Eye Massager" on the mobile phone to connect. The indicator light will be solid blue after connection is successful.

- When connected, the voice broadcasts "Bluetooth Connected".
- When disconnected, the voice broadcasts "Bluetooth Disconnected".

FR

IT

ES

Long press " " for 3s to enter/exit pairing mode.

## FAQ and Troubleshooting

Problem	Possible causes	Actions
There is a sound mechanical or air sound coming from device	Normal sound of air pump inflating, motor working, or air pressure	No need to do anything
Device turns off automatically. Device does not turn on. Buttons are not responsive. Air pressure is weak.	1.Battery is Low. 2.The 15 mintues timer.	1.Charge the device. 2.Turn on the massager again.

If there are other problems, turn off the device and contact us at [support@renpho.com](mailto:support@renpho.com). RENPHO is not responsible for any damage or injuries that occur if you attempt to disassemble or repair the product yourself.

## Cleaning and Maintenance

Attention

Please turn off the power when cleaning, do not use corrosive detergents such as benzene or diluents to clean the product.

Do not expose the device to excessive sun or hot temperature.

## Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

**<https://renpho.com/pages/warranty-terms-and-conditions>**

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

## Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

**TEL: +1(844) 417 0149 (US ONLY)**

Monday-Friday 9:00AM-4:30PM

**1800 MY RENPHO (1800 69 7367) (AU ONLY)**

Monday-Friday 9:00AM-5:00PM (GMT+10)

**+44(0)7434-666088 (UK ONLY)**

Monday-Friday 9:00AM-5:00PM(GMT)

**Email: [support@renpho.com](mailto:support@renpho.com) (US&CA)**

**[support-au@renpho.com](mailto:support-au@renpho.com) (AU)**

**[support-eu@renpho.com](mailto:support-eu@renpho.com)(UK&EU)**

**\*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.**

### Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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## FCC Regulatory Compliance

---

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

### RF Exposure Compliance

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

### ISED Regulatory compliance

This device contains licence-exempt transmitters that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s).

Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

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