

RENPHO

RF-FM062

FOOT MASSAGER

USER MANUAL




Introduction

Thank you for choosing the RENPHO foot massager with kneading and heating function.

Please read the user manual carefully before use and keep it in a safe place for future reference. If you need assistance, please contact our support team.

Product Diagram:

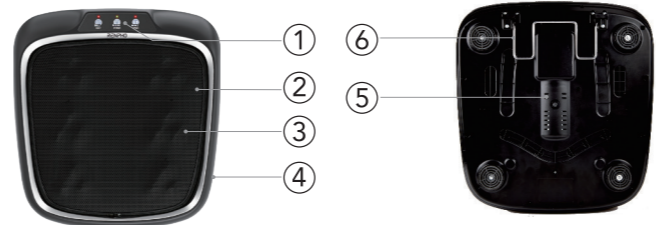

- | | |
|--|-----------------------|
| 1. Control panel buttons: Heat, Power, Speed | 4. Power interface |
| 2. Massage nodes | 5. Heat emission hole |
| 3. Washable cover | 6. Adjustable stand |

Function

HEAT: Press to turn on/off the heat function

POWER: Press to turn on/off the massager

SPEED: 2 speeds

Removable mesh: Unzip to remove and clean; machine-washable

Adjustable stand: Adjusts the height and angle of the foot massager

-01-

Product Operation Guide

1. Connect the plug to a power source.
2. When the power indicator is red, the massager is on standby. When it is green, the massager is working.
3. Press the "POWER" button to activate the massager. The first and the third massage nodes on the left foot will rotate counter-clockwise, while the second massage node rotates clockwise. The respective nodes rotate in the opposite direction for the right foot.
4. Press the "HEAT" button to activate the heat function. Press it again to turn off the heat function.
5. Press the "SPEED" button to adjust the rotation speed. Flashing red light indicates slow speed. Solid red light indicates fast speed.
6. The foot massager will automatically turn off after 15 minutes to avoid over-heating.
7. To remove the cover, simply unzip. Cover is machine-washable.
8. As a safety feature, the motor will auto shut off when the temperature gets too high, either through operating for too long or having too heavy of a load. When the temperature cools, the motor will operate normally.

Note: We recommend you limit usage to 30 minutes. Longer massage time may cause mild pain.

WARNING

1. This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and carefully aware of the risks of use.
2. Children must not play with the device.
3. Cleaning and maintenance must not be performed by children unless supervised.
4. If the mains connection cable of this device is damaged, it must

-02-

FCC Statement

This device complies with ICES-003 Rules and US Electronic Code of Federal Regulations 47 CFR FCC Part 15; Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Statement NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener. DO NOT hang the unit by the power cord.

US Importer: JOICOM CORPORATION

14129 The Merge Street, Building 3 Unit A, Eastvale, CA 92880

Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.

Qianhai Complex A201, Qianwan Road 1, Qianhai Shenzhen-Hong Kong Cooperation Zone, Shenzhen, P.R.China

Made in China
NA3

-05-

be disposed of. If it cannot be removed, the device must be disposed of.

5. Do not insert any objects in the device apertures or in the moving parts. Ensure that the moving parts can move freely at all times.
6. Do not force or jam any part of your body, such as your fingers, between the moving massage elements or their holders.
7. Only to be used with the power supply unit provided with the appliance.
8. Do not stand on or in appliance. Use only while seated.

CAUTION

- For environmental reasons, do not dispose of the device in the trash.
- Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with CE Directive- WEEE (Waste Electrical and Electronic Equipment).
- If you have any questions, please contact the local authorities responsible for waste disposal.

Troubleshooting Guide

Common problems	Cause	Solutions
The device does not work	1. Disconnected power or loose plug causing poor contact. 2. The motor auto shuts-off due to the device being on too long activating	1. Re-plug into the power source. 2. Use the massager after the motor cools off naturally.
The heat function is not working or not as expected	1. The heat function is not on. 2. Can't feel the heat.	1. Press the "HEAT" button. The red indicator light should turn on. 2. It takes about 7-10 minutes to fully heat up. Please be patient as this is to avoid the temperature rising too fast.

-03-

Package included

- 1 x RENPHO Foot massager
- 1 x User manual
- 1 x Power Adapter

Product Specifications

Model: RF-FM062
Rated Voltage: 12V $\overline{=}$ 2.5A
Rated Power: 30W
Default Timer: 15 mins

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:
<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

☎ TEL: +1(844) 417 0149 (US ONLY)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-5:00PM (GMT+10)
✉ Email: support@renpho.com (US&CA)
support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

-04-